

Image Quantification Toolbox Installation Info

Fiji Plugins:

- Java-8 (enabled by default)
- ImageJ-MATLAB plugin (vital)
- OrientationJ plugin (needed for fiber orientation)
- ClearVolume plugin (needed for built-in 3D viewer)

MATLAB Toolboxes:

- Simulink (needed for cell analysis)
- MATLAB Computer Vision Toolbox (needed for cell analysis)

Installing Fiji Plugins:

Step 1:

Open the installation of Fiji you would like to use inside the toolbox.

Step 2:

Navigate to Help > Update...

Step 3:

Click "Manage update sites".

Step 4:

Select the Fiji Plugins listed above and click "Close".

Installing MATLAB Toolboxes:

Step 1:

Open MATLAB

Step 2:

Navigate to Apps > Get More Apps

Step 3:

Search for the toolboxes listed above and add them to MATLAB.

Installing the Image Quantification Toolbox:

Step 1:

Install the dependencies above.

Step 2:

Launch InstallerApp.mlapp.

Step 3:

Select your desired installation option and click “Install”.

Step 4:


Navigate to the folder containing the version of Fiji you would like to use and click “open” in the file explorer.

Step 5 (optional):

If you selected “Custom” for your installation location navigate to the folder you would like to house the toolbox.

Step 6:

Click “Browse for folder” button above “Current Folder” and navigate to the folder containing the toolbox.

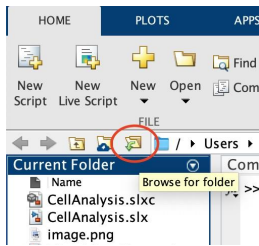

Step 7:

Please launch ImporterApp.mlapp to begin using the toolbox.